

## US300

### Ultrasonic Cleaning for Cannulated and Robotic Instruments



#### Models

US 300                      W x D x H   35.1" x 26.6" x 41.2" (lid closed)

#### Standard Features:


- LCD Color Display
- Cleaning tank in AISI 316L stainless steel
- Three water connections CW, WW, DI
- Dosing pump for detergents
- Flow meter monitoring precise detergent dispensing
- Detergent level sensor
- Automatic filling and draining of chamber
- IR-proximity sensor for automatic opening / closing of lid
- RS232 interface for external printer
- USB port for cycle data download
- Built-in heaters for thermal rinse at 195°F
- Drain discharge cool down
- Two rotating spray arms for instrument surface cleaning
- Powerful ultrasonic generator and transducers
- Steelco Active Flushing Technology with 20 flush ports for cannulated instruments
- Integrated detergent storage compartment
- Two US300503R trays with adjustable and removable instrument supports
- Forty instrument connection sets
- Compatible with Intuitive DaVinci Xi-/Si Robotic Instruments
- Mobile with four casters
- Thermal self-disinfection cycle



#### **STEELCO USA Inc.**

7354 Central Industrial Drive West Palm Beach FL 33404 USA

Ph. 561 791 8313   Fax 561 660 5382   www.steelcogroup.com   e-mail: contact@steelco-usa.com

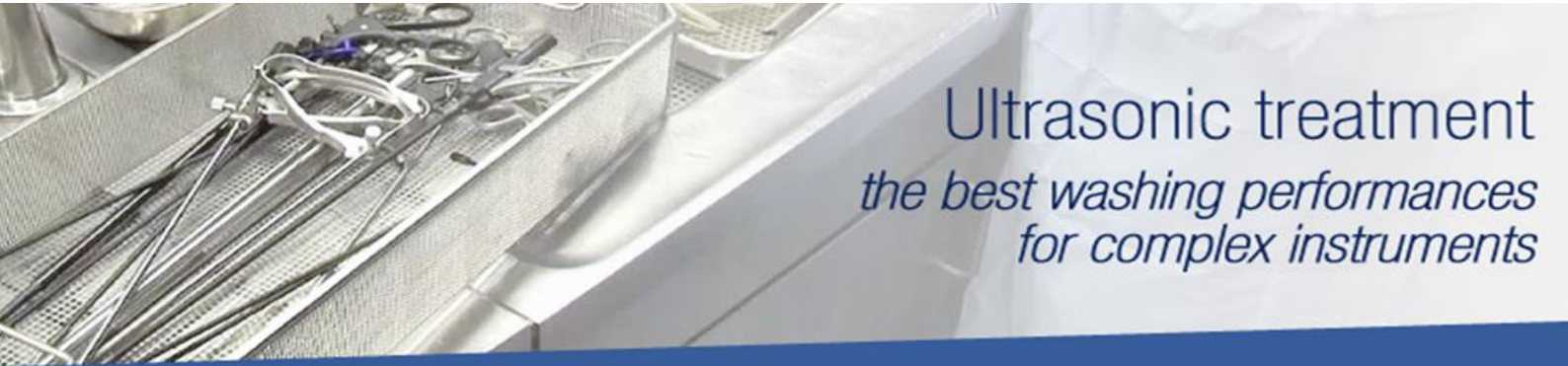

Ultrasonic treatment  
*the best washing performances  
for complex instruments*

### **Available Configurations:**

- Electrical connection 208V, 3-phase


### **Cycle Times:**

#1 Ultra Rapid	28 min including cleaning and rinsing
#2 Standard	43 min including cleaning, rinsing and thermal disinfection
#3 RoboPro	65 min including pre-cleaning, cleaning, rinsing and thermal rinse

### **Options:**

- ST2 Integrated Printer

### **Instrument Cleaning Inserts:**

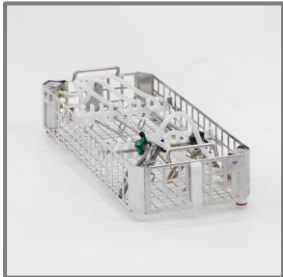

#### **US300503R Instrument Tray**

- Twenty flush ports for cannulated instruments
- Adjustable and removable instrument supports
- Accepts instruments up to 23.5" total length

#### **US300518R Connection Kit DaVinci Si-Series**

#### **US300530R Connection Kit DaVinci Xi-Series**

### **Consumables:**



#### **US405033 MediClean forte 5L**

- Liquid concentrate detergent  
Excellent material compatibility, universally applicable. Residues of dried and denatured blood are reliably removed. A redeposition of protein residues is prevented.